

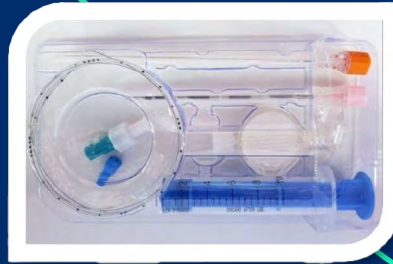


Epidural Kits

by SURU®

**Mini Pack System I**

Item Code # 3038.1

**Combined Spinal Epidural Kit System II**

Item Code # 3038.2

**NR-FIT Epidural Kit I**

Item Code # 3038N.1

**NR-FIT Combined Spinal Epidural Kit II**

Item Code # 3038N.2

SURU® Epidural Kits

Product Features

- Epidural catheter kit are designed for short term & long term continuous anaesthesia & pain therapy
- Sterile, individually packed in blister packing
- Also available NR-Fit



Australia
1800 358 101
sales.au@westlab.com
www.westlab.com.au

Canada
1 877 822 1455
sales.ca@westlab.com
www.westlab.com

New Zealand
+64 9 553 3677
sales.nz@westlab.com
www.westlab.co.nz

United States
646 568 5391
sales.us@westlab.com
www.westlab.com

Epidural catheter

- Epidural catheter available with closed soft tip or closed normal tip, with three lateral eyes
- With radio-opaque line
- Length : 1000 mm

Epidural / tuohy needle

- Tuohy type / epidural needle is provided with clear depth marking for accurate insertion depth reading
- Length of epidural needle : 80 mm
- Polished inner bevel edge minimize the risk of catheter shearing
- Size : NG 16, 17, 18, 19

Lor syringe

- Specially designed plastic syringe is provided for clear identification of epidural space
- LOR syringe is suitable for air and saline technique
- Syringe capacity: 10 ml

Threading assist guide

- Fit securely into the needle and assists in easy advancement of catheter

Detachable connector

- Separate twist lock connector is provided for safe and secure attachment to the catheter for convenience of the procedure

Available Configuration:

Epidural Anaesthesia Kit Consist of :

Epidural Kit

System I Kit Consist of :

- Epidural Catheter
- Epidural Needle
- LOR Syringe
- Hydrophilic Filter
- Detachable Connector
- Thread Assist Guide

Epidural Spinal Combined Kit

System II Kit Consist of :

- Epidural Catheter
- Epidural Needle
- LOR Syringe
- Hydrophilic Filter
- L P Spinal Needle
- Detachable Connector
- Thread Assist Guide

For more information:

<http://www.suru.com/Anaesthesia-Respiratory-Care/Epidural-Kit>

Australia

1800 358 101
sales.au@westlab.com
www.westlab.com.au

Canada

1 877 822 1455
sales.ca@westlab.com
www.westlab.com

New Zealand

+64 9 553 3677
sales.nz@westlab.com
www.westlab.co.nz

United States

646 568 5391
sales.us@westlab.com
www.westlab.com